
USER'S MANUAL

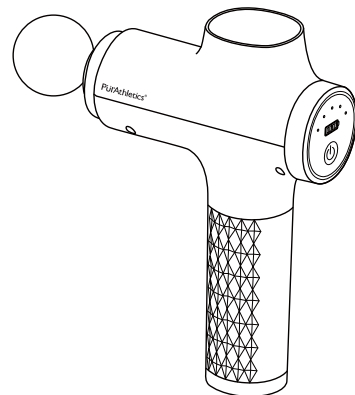Product Specification

Percussion Pulse Massager

Size : 168.5 mm X 238.5 mm X 63 mm

Net weight : ≈ 680 g

Battery capacity : 1800mAh X2

Charging interface : USB type-C

Charging time : ≈ 3.5 hours



Maximum (blocking) current : about 7.0A

Gear/speed: 6 Gear/1800 ~ 3000 RPM

Stroke: 6mm


Maximum Power: about 60W

Maximum operating noise: about 50dB

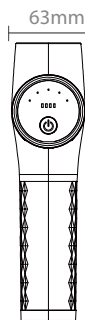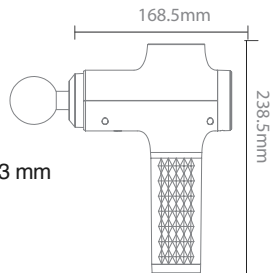
Parts description and functionality

- ① Quick release massage attachment
- ② Grip
- ③ Gear indicator
- ④ Power indicator
- ⑤ On/off switch operating mode
- ⑥ USB type-C charging connector

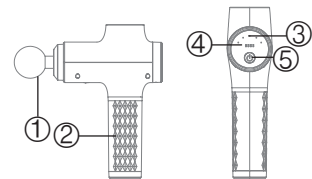Accessories

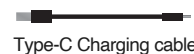Quick guide usage and functionality

1. Massage head installation.

Turn the massager off before installing the massage head. Choose the suitable massage head according to the area you need to massage. Carefully insert the massage head into the receiver.

2. Long press the on/off button for about 1.5 seconds, the mode indicator and the power indicator will light up. Short press the button to choose modes. By short pressing the power button, mode 1 will turn on and the rest of the indicator lights continuously flash.

This device is equipped with smart mode functionality. Pressing the massage head firmly enables the head to massage deeply with alternating speed pattern.

Long press the power button to turn it off.

3. This device has a default of automatically stopping every 5 to 7 minutes. To continue your massage short press the power button to extend your session to another 10 minutes.

4. This state of the art device is equipped with a short press button that shifts to various massage modes from low to high cycle. The following are the suggested mode usage descriptions:

- Mode 1 suitable for warming up muscles
- Mode 2 suitable for relaxing the fascicles
- Mode 3 suitable for lactic acid removal
- Mode 4 suitable for deep massage
- Mode 5 suitable for professional massage
- Mode 6 all the 5 mode lights flashing simultaneously

This is an automatic frequency conversion, depending on the pressure, it automatically adjusts the strength and speed of massage.

5. Charging instructions:

- (1) When the power is insufficient, the first power indicator lamp flashes.
- (2) The corresponding power indicator light flashes when the unit is being charged.
- (3) Once the charge is complete, all four-unit power indicators will turn on.

Important safety tips

1. This device is suitable for 5V/2A chargers or partial fast charging adapters.
2. Fully charge the battery before first use.
3. Please make sure the input voltage is the same as the adapter voltage of the device before charging.
4. Do not change the massage attachments while the device is running.
5. Do not immerse the product in water.
6. Do not apply vibration to the head area, groin area, chest or areas of the body where surgery was performed.
7. Do not use the device as an alternative to treatment or adjuvant therapy. If you have pain or are currently undergoing treatment or using other medical devices, consult a specialist before using this device.
8. In case of malfunctions, contact a technician, do not disassemble the product yourself.

WARNING

In the following cases, as well as for any health problems, consult an expert before using the product.


1. Pregnancy, diabetes, and complications such as nerve and retinal damage.
2. The presence of a pacemaker, recent surgery, epilepsy or migraine, prolapse of the lumbar inter vertebral disc, displacement of the vertebrae.
3. Recently underwent joint replacement surgeries, IUD placement, metal implants, or other health concerns.

ENVIRONMENTAL SAFETY

This product contains a lithium battery and must not be disposed of with normal household waste and must be disposed of safely at a designated hazardous waste or electronic product collection point.

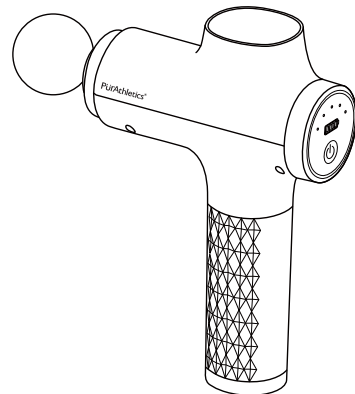By using this percussion massager, you acknowledge and accept the risks associated with its use. The company, distributors, and resellers shall not be held liable for any injuries, damages, or health issues resulting from the improper use or misuse of this product. If you experience any discomfort or adverse reactions while using the device, discontinue use immediately and seek medical attention.

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