PürAthletics®

Toll Free Consumer Support • Numéro de téléphone sans frais 1-877-415-3425  
Zenation Athletics, Zensation Athletics Rejuvenate, PürAthletics, PürEarth and PawsAthletics  
products are distributed by TriMax Sports, Vancouver B.C, Canada V5X 4V6  
Les produits Zensation Athletics, Zensation Athletics Rejuvenate, PürAthletics, PürEarth et  
PawsAthletics sont distribués par TriMax Sports, Vancouver, C.-B. Canada V5X 4V6  
email/courriel: service@zenathletics.com | 9 am to 5 pm (PST)  
Mon - Fri | 9 à 17 h (HNP) Lun - ven

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Product After-Sales Service

Thank you for choosing this product. In order to provide you with better after-sales service, we kindly request your cooperation as follows:

1. When purchasing, please ask the seller to issue a sales certificate as the basis for warranty;
2. Please fill out this card in detail;
3. Please carefully read the maintenance regulations and each clause;
4. If you need service, please call the customer service hotline directly;

Maintenance Regulations:

First, during the warranty period, this product comes with a one-year free warranty for the entire machine under normal use.  
(From the date of purchase, the purchase date shall prevail).  
Second, during the warranty period, if one of the following situations occurs, material costs or service fees will be charged at your discretion, We apologize for any inconvenience caused!  
1. Consumers who cause product damage due to improper use, maintenance, or storage;  
2. Damage caused by disassembly by non warranty repairers;  
3. Appearance that naturally stains due to use;  
4. Damaged for commercial use;  
5. Damage caused by force majeure is not covered by the warranty.  
Third, this card cannot be reissued, please keep it safe.  
Fourth, for products that do not fall within the scope of repair or exceed the warranty period, we still warmly serve you and only charge material cost fees.

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Certificate Of Conformity  
Certificat de Conformité

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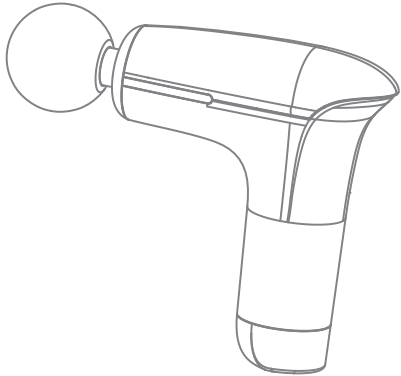
USER'S MANUAL | MANUEL D'UTILISATION

Dear User

Thank you for choosing to use this product. This product is a well-designed and high-quality household appliance product. For your convenience, please read the manual carefully and operate according to the instructions in this book. From beginning to end, our company will provide you with high-quality service. No matter what problems you have during use, please feel free to call our company's service hotline or contact our local service center. Thank you again for choosing to use our company's product. Due to product improvements, the product you received may not be completely consistent with the illustration in the manual. We apologize for any inconvenience caused.

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Matters Needing Attention

The following precautions must be followed during the use of this product:

- Do not use this product while charging.
- Do not use this product on the injured area.
- Never overload the product during operation.
- Do not disassemble the product by yourself.
- Do not attempt to modify the product yourself.
- When removing the massage head, do not insert your fingers into the connection of the massage head.
- When the product is running, please do not use the product to point at people.
- Do not massage one area for too long to avoid discomfort.
- Do not soak or place the product in a damp environment.
- Do not use this product to massage parts of the body that are prone to injury.
- Cancer patients or patients with major illnesses are not allowed to use this product.
- Please do not dispose of it casually or use it in an unreasonable way.
- Do not operate this product without the user's attention.
- Do not wear loose clothing when using this product.
- Do not use this product to massage the head.
- Do not use it outdoors.
- Pregnant women, implantation of pacemakers, fractures, osteoporosis, implantation of any metal medical devices, etc., should not use this product.
- After 10 minutes of use, please shut down for at least 5 minutes.
- Please ensure that the ventilation holes are unobstructed.
- Please keep your hair and fingers away from the piston rod and ventilation hole to avoid getting caught.
- If you feel any problems with your body or skin, please consult your doctor and follow their instructions.
- To avoid excessive stimulation of muscles and nerves, it is recommended that the massage time should not exceed 10 minutes at a time.
- This handheld fascia gun is equipped with a heating component, and should be used by people who are not sensitive to heat to avoid burns;
- The use of unauthorized or incompatible power sources, chargers, or batteries may cause fire, explosion, or other hazards;
- The product is a non-medical device and cannot be used for medical purposes. If it is intended to be used for the recovery of muscle injuries, it is recommended to consult a doctor.
- Children are not allowed to play with equipment.
- Children are not allowed to perform cleaning and user maintenance without supervision.

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**Warning 1:** When using this product, users should use appropriate force and avoid forceful pressing that may cause the device to stop (overload). Excessive "overload protection" can cause destructive damage to the device.

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Cleaning And Maintenance

- Use a soft cloth to wipe the surface of the product.
- It is prohibited to immerse the entire product in water or other liquids during cleaning.
- Use corrosive cleaning agents such as benzene or diluents to clean the massage gun.
- Sweep away dust. Keep the massage gun clean and away from damp places.
- Store the massage gun in a cool and dry place out of reach of children.

How to remove the built-in rechargeable battery  
This product contains a rechargeable battery. When consumers decide to dispose of this product, please make sure to remove the battery from the product:  
Before disposing of the appliance, the battery must be removed from the appliance;  
When removing the battery, the appliance must be powered off;  
Batteries should be disposed of safely.

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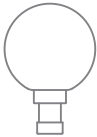
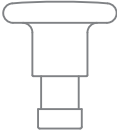
Product Usage Method

- Press and hold the power button for 3 seconds to start the product in a low-speed state.
- By analogy, each press can achieve a suitable speed and reach a massage frequency that suits oneself. This adjustment will repeat the cycle.
- During operation, press and hold the power button for 3 seconds, and the machine will stop working. All lights will turn off.
- When the product is being charged, the machine cannot work, and the light strip will flash with a blue blinking light. When fully charged, the blue light will stay on.
- The machine will automatically stop working continuously for 10 minutes (this setting is based on protecting the motor to extend the product's service life). After stopping, please restart it for 3-5 minutes.
- When the user uses too much force to exceed the product's maximum load, the machine will stop rotating. After the state lasts for 5-10 seconds, the machine will automatically shut down.

Attention

- This product can be connected to a laptop, phone adapter, power supply, or any power supply that provides 5V DC through a Type-C data cable Line charging.
- The charging cycle varies depending on the power supply.
- After charging, please remove the power cord from the power supply

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