
SETTING UP YOUR **CORESENSE SMART SCALE** INSTRUCTIONS

OVERVIEW
Leaone stands as your ultimate intelligent hardware management platform, offering comprehensive features and seamless connectivity. With the user-friendly Lelife App, effortlessly establish a quick and convenient connection between your mobile device and the CoreSense Smart Scale. Conduct multiple tests, securely store data, and enjoy unparalleled convenience, all at your fingertips.


DOWNLOAD AND INSTALL THE APP FOR FREE
Get started in seconds! Scan the QR code or search 'Leaone' in your app store. Enjoy hassle-free installation and explore the full capabilities of our platform.

REGISTRATION
1. Setting up the account: choose how you want set your account using the following options, mobile phone or email. The information you provide will be used later when you log in.
2. Account Verification: Click the "Send" button to dispatch the verification code to the registered account, ensuring its accuracy. The transmission time for the verification code is set at 120 seconds. If the code is not received within this time frame, simply click "Send" again.
3. Password Setup: To uphold account security standards, passwords must adhere to the following criteria: 6-16 characters in length, including numbers and letters.
4. Click on register to complete the process.

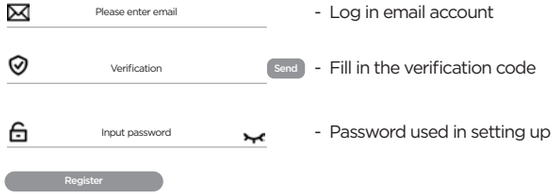- Log in email account
- Fill in the verification code
- Password used in setting up

TO LOG IN
Enter your account + password.
You can also login using third party apps including wechat, facebook and in the app store using your apple ID.

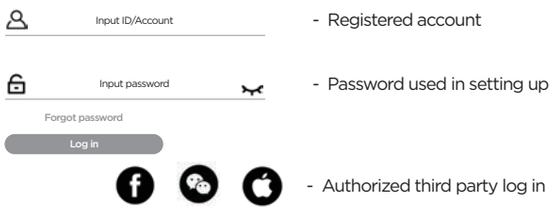
- Registered account
- Password used in setting up
- Authorized third party log in

***Registration entails the user's agreement to abide by the Service Agreement and Privacy Policy.**


RESETTING LOST PASSWORD
You can easily reset your forgotten password by following these steps:
1. Enter the account that you want to retrieve.
2. Send the verification code. Click the SEND button to send the verification code to the input account number and verify the account number. The transmission time of the verification code is 120 secs. If you don't receive the verification code, please click Send again.
3. Reset the password using the same steps when you registered.
4. Once the new password has been completed, click Submit Reset to log in to the app using the new password.

CREATING A NEW PASSWORD

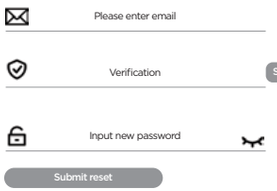
- Log in email account
- Fill in the verification code
- Type in the new password

ACCOUNT INFORMATION
To accommodate the range of supported devices, users are required to configure an avatar, nickname, birthday, and gender affiliation. This additional information is utilized for data calculations in certain devices and for analyzing measurement data. Users can adjust settings such as nickname, birthday, gender, and height, alongside electrode measurement settings, and then proceed to calibrate the equipment accordingly.

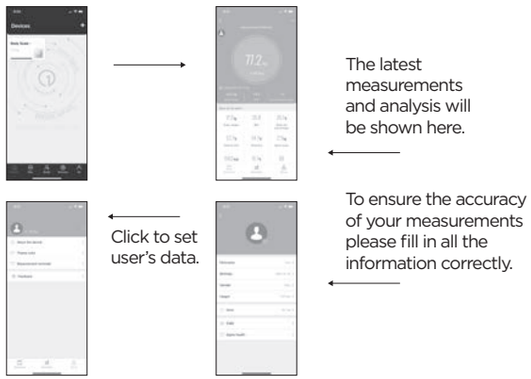PAIRING/CONNECTING TO THE DEVICES
There are 2 ways to pair the equipment and the devices either choose nearby equipment or add them manually.
1. In nearby equipment, open the nearby equipment page to automatically search the matching equipment close by, click the equipment found to pair.
2. In manual addition, select the equipment to be added on the list and add according to the operation guidelines.

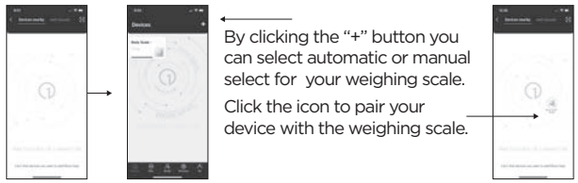DATA FUNCTION
The data function serves as a comprehensive display of user-related measurement data summaries. For instance, if a user utilizes two different body fat scales, both producing weight data for the same individual, this information is consolidated on a single user profile within the data page. Users can effortlessly click on the corresponding data type to access a detailed summary of all recorded data points.

FRIEND FUNCTION
Users can easily compare data between themselves by connecting with friends.

HOW TO ADD FRIENDS
There are two ways to add your friends, search to add the account ID/account or scan their QR code.
Method 1: Search for your friend's ID/account: open the "Search ID/account" function to enter the ID or account of your friends who need to be added, and search to add.
Method 2: Scan add, open the scan function to scan the QR code of your friend you want added.

VIEWING YOUR FRIEND'S DATA
After searching or scanning, the user profile page will open. Navigate to "Data Open User" within "Privacy Protection," where you can select which family members you'd like to share with friends. Only the chosen family members will be visible to your friends. If no family members are selected, friends will not have access.

FAMILY MANAGEMENT
In LEAONE, numerous devices support multiple family members' usage. With the ability to switch between different families, you can conveniently store data under separate family profiles.

CREATE/MANAGE FAMILY MEMBERS
To do this, open MINE and select My Family to manage and create the family . Choose the corresponding family to edit the basic information. Choosing Add Family creates a new family and therefore it needs setting up a new profile used for the calculation, analysis and data gathering.

Note: The relationship is a unique identification between family and the main account holder and cannot be modified once the information has been inputted.

REMOVE A FAMILY MEMBER
By selecting the family member, you can click on the trash bin to delete the account. Deleted information will no longer be retrieved.

SETTINGS
1. Click settings to go to Display Mode. You can choose your preference of display to Grid or List.
2. Click the language menu to change the language you prefer. The default is auto, English.

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